



WIRELESS CONNECTIVITY

resoundpro.com

Extend your
patients' hearing
experience

ReSound



rediscover hearing

Empower your patients with all the benefits of Smart Hearing™

ReSound Smart Hearing adapts to the way people want to live today. It works so intuitively that your patients will find surprising new ways to use their hearing aids to enrich their lives. With an entire ecosystem of wireless connectivity for streaming sound, on-the-spot fine-tuning and personalization, your patients will hear more, even in noisy places, so they have more confidence to join in.

Whether it's a phone call, a movie or a quiet conversation in a busy restaurant, every sound is adjusted to each individual's hearing loss to deliver the most natural, effortless hearing experience possible.

Smart Hearing. A whole new world of wireless connectivity that brings your patients closer to what they love.



Hearing aids that work like wireless stereo headphones

Our 2.4 GHz technology makes it possible to connect ReSound Smart Hearing aids directly to an iPhone®, iPad® or iPod® touch.

This means your patients can stream music, make phone or video calls, listen to podcasts and audiobooks, and even get turn-by-turn directions around an unfamiliar city.

If patients don't have an iPhone, they can still stream stereo sound for calls and music from their Android smartphone or any other Bluetooth-enabled phone with the ReSound Unite™ Phone Clip+.

Made for iPhone connectivity. The simplest and most direct way to make hearing aids work like wireless stereo headphones.



Made for
iPod iPhone iPad



WHY RESOUND WIRELESS TECHNOLOGY LEADS THE INDUSTRY

ReSound is the only manufacturer with a comprehensive range of seamlessly integrated 2.4 GHz wireless accessories.

For your patients, this means excellent sound quality and robust connections, all with the same pairing protocol and low latency.

So your patients don't have to choose between the best hearing aid, the best wireless accessories or Made for iPhone connectivity.

They can have it all.

Extend your patients' hearing and access to speech

The ReSound Micro Mic and ReSound Multi Mic are discreet, clip-on microphones that increase signal-to-noise ratio so your patients can understand speech in challenging listening environments.

RESOUND MICRO MIC

Extremely easy to use, the ReSound Micro Mic pairs to any ReSound Smart Hearing aid* in seconds. When your patient clips it on to the clothing of the person they want to hear, they can enjoy conversation even in background noise. With a range of over 80 feet, the ReSound Micro Mic is ideal for use on shopping trips, in busy restaurants or for talking in the car.

RESOUND MULTI MIC

The ReSound Multi Mic can do everything the ReSound Micro Mic can and more. Placed flat on the table, the ReSound Multi Mic automatically turns into a table microphone so your patients can hear what everyone around the ReSound Multi Mic is saying. The ReSound Multi Mic also connects to FM and loop systems**, and has a mini-jack input for streaming audio from virtually any device with a headphones output.



ReSound Micro Mic significantly improves signal-to-noise ratio in noisy situations, delivers exceptional speech understanding and extends your patient's hearing over 80 feet*.

ReSound Multi Mic, with the same hearing range over 80 feet, doubles as a table microphone, connects with loop and FM systems**, and has a mini-jack input for streaming audio from a computer or mobile device.

* Compatible with all ReSound wireless hearing aids
** FM receiver required



ReSound Multi Mic

Place it on the table and it automatically switches into an omni mode, optimized for picking up sounds around the ReSound Multi Mic.

Plug into virtually any mobile device or computer and stream sound.

Connect to an FM system** and stream sound.

Rechargeable 11 hour battery life.

Switch between microphone, mini-jack input, FM and loop systems at the push of a button.

Robust clothing clip rotates 360° for optimal voice pick-up.





Home just got a lot more entertaining

The ReSound Unite™ TV Streamer 2 makes it easy for the whole family to watch TV together at a volume that's comfortable for everyone.

EASY TO CONNECT

Your patients can simply plug the ReSound Unite TV Streamer 2 into their TV, music system or computer. Once it is paired to their ReSound Smart Hearing aids*, they can enjoy crystal clear stereo sound, and listen at their volume while others around them listen at theirs. If your patients leave the listening range and come back within five minutes, their hearing aids will automatically reconnect.

WATCH AND CHAT AT THE SAME TIME

Your patients can use their ReSound Smart™ app or ReSound Unite Remote Control 2 to adjust the volume of their hearing aid microphones and the TV independently. This means they can chat with whoever is sitting next to them, or they can mute their hearing aid microphones altogether and concentrate on the TV program.



ReSound Unite™ TV Streamer 2

Your patients can stream stereo sound from their TV, computer or music system directly to their hearing aids.



* Compatible with all ReSound wireless hearing aids.



Connect your patients with crystal clear calls

The ReSound Unite™ Phone Clip+ is a discreet device that connects wirelessly to Bluetooth-enabled phones. It streams the caller's voice directly to your patient's hearing aids, while the microphone on the ReSound Unite Phone Clip+ picks up their own voice. If background noise begins to intrude, there's a mute button that deactivates the hearing aids' microphones.

STREAM AUDIO

The ReSound Unite Phone Clip+ is not just for phone calls. Your patients can stream any audio coming from their phone, such as music. With the ReSound Unite Phone Clip+, their hearing aids work like wireless stereo headphones.

CONTROL

They can even use the ReSound Unite Phone Clip+ as a basic remote control to adjust hearing aid volume and change programs.

PUT CONTROL AT THEIR FINGERTIPS

For patients who find it tricky to use the buttons on their hearing aids, the ReSound Unite™ Remote Control 2 gives them intuitive control over volume and program settings. Its large buttons also give them access to all of their ReSound streaming accessories, and they can check their settings and connections on the backlit display.



ReSound Unite™ Phone Clip+


Patients can make crystal clear phone calls, mute background noise while they talk, or stream music, podcasts or any other audio from their phone.



ReSound Unite™ Remote Control 2

Patients can adjust the volume or mute their hearing aids, switch programs, and see all settings at a glance.

Award winning app for a personalized hearing experience

 With the award-winning ReSound Smart™ app, your patients can personalize their hearing experience from an iPhone or Android smartphone, or even an Apple Watch.

DISCREET ADJUSTMENTS

With the ReSound Smart™ app, they can adjust volume, treble and bass wherever they are. Use the Sound Enhancer features for better speech understanding, spatial awareness and listening comfort. Your patient can create favorite programs, geotag them, and have them activated automatically as they enter a new location, all without having to touch – and draw attention to – their hearing aids.

PERSONAL WIRELESS NETWORK

Your patients can control their other ReSound streaming accessories, too. They can turn up the volume of their ReSound Micro Mic or ReSound Multi Mic when they want to focus on what the person wearing it has to say – and turn it down when they don't. And, when their favorite program comes on TV, they can switch to their ReSound Unite TV Streamer 2 with a couple of taps on their smartphone.




Enhance sound in challenging sound environments

The ReSound Smart app includes three advanced features in the Sound Enhancer.



SPEECH FOCUS
This enables your patients to make three directional beamwidth adjustments.

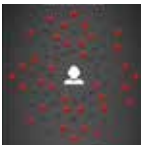


COMFORT IN NOISE
Your patients can reduce noise levels in loud environments so listening is more comfortable and understanding speech takes less effort.




COMFORT IN WIND
Your patients can adjust the wind filter to reduce the amount of intrusive wind noise for maximum listening comfort.



TINNITUS MANAGEMENT
If you activate Tinnitus Management in your patient's hearing aids, they can fine-tune the sound therapy with adjustments to pitch and modulation, or they can choose between six different nature-inspired sounds.

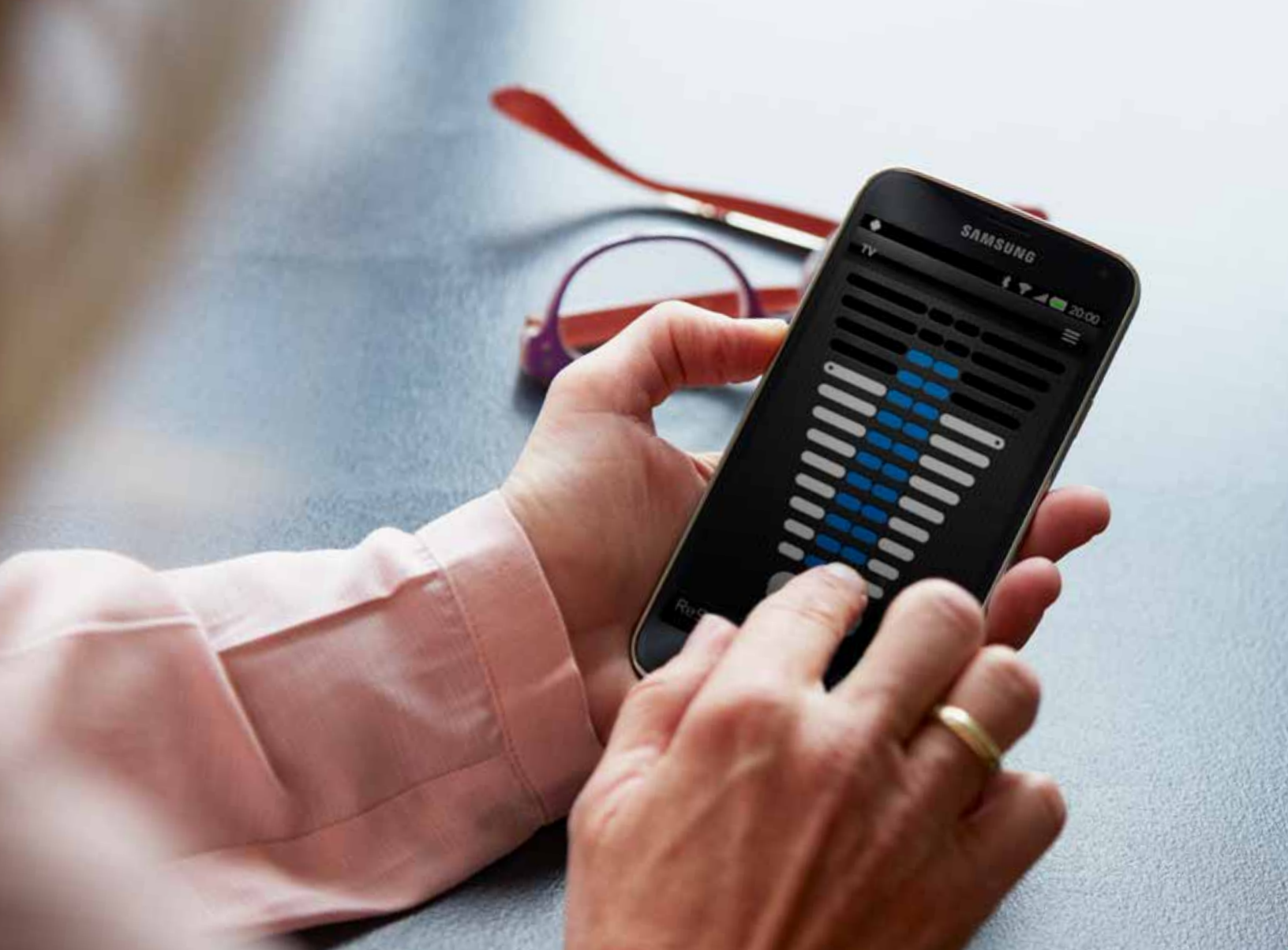

Discreet control for all wireless hearing aids

 With the ReSound Control™ app, adjusting volume or changing programs is as easy as a tap and a swipe.

EASY SETUP WITH PHONE CLIP+

If your patient has a wireless hearing aid or smartphone not compatible with the ReSound Smart™ app, the ReSound Control app is a good and discreet solution for adjusting volume, changing programs and controlling streaming accessories such as the ReSound Micro Mic, ReSound Multi Mic or ReSound Unite TV Streamer 2.

The ReSound Unite Phone Clip+ is needed to use the ReSound Control app.



ReSound Control app. Your patients can change volume, switch programs and balance the volume of their wireless hearing aids with their ReSound streaming accessories.

Easy overview

WITH RESOUND SMART HEARING AIDS, YOUR PATIENT GETS ACCESS TO:

- An iPhone, iPad or iPod touch
- A Phone Clip+
- Any combination of three streaming accessories
- A Remote Control 2





PHONE AND APP OPTIONS



Phone	Smart Hearing aids	Wireless hearing aids
IPHONE (IOS7 OR LATER)	Phone call	Made for iPhone
	App use	ReSound Smart
ANDROID PHONES (SELECTED MODELS)*	Phone call	Phone Clip+
	App use	ReSound Smart
OLDER IPHONE OR OTHER ANDROID* PHONES	Phone call	Phone Clip+
	App use	ReSound Control
OTHER BLUETOOTH-ENABLED PHONES	Phone call	Phone Clip+
	App use	Phone Clip+

NOTE: Other combinations are also possible
* See www.resound.com/compatibility for list of phone models compatible with ReSound Smart app

CONTROL OPTIONS



	ReSound Smart	ReSound Control*	Phone Clip+	Remote Control 2
VOLUME AND PROGRAM FUNCTIONS				
Volume adjustments	Mono-/binaural	Mono-/binaural	Binaural	Mono-/binaural
Hearing aid program selection	Direct select	Direct select	Toggle (1-2-3-1...)	Toggle (1-2-3-1...)
Streamer selection	Direct select	Direct select		Toggle (1-2-3-1...)
Balance hearing aid volume vs. streaming volume	●	●		●
Hearing aid mute	●	●	●	●
OTHER FUNCTIONS				
Phone operation			●	
Program/streamer name customization by user	●	●		
Bass/treble adjustment	●			
Display	●	●		●
Create favorite programs with or without geotags	●			
Sound Enhancer adjustments**	●			
Tinnitus Management adjustments	●			
Find your lost hearing aid	●			

* ReSound Unite Phone Clip+ is needed to use the ReSound Control app
** Available in ReSound LINX² 9 and ReSound ENZO² 9

STREAMING OPTIONS



	Phone Clip+	Multi Mic	Micro Mic	TV Streamer 2*
STREAMER FUNCTIONALITY				
Power Source	Rechargeable – portable	Rechargeable – portable	Rechargeable – portable	Wall plug-in A/C adapter – stationary
Audio Quality	Stereo (A2DP)	Mono	Mono	Stereo
Voice Pickup	Directional microphone	Directional/Omnidirectional	Directional microphone	
Telecoil		●		
Line-in		●		
FM compatibility (FM receiver required)		●		

* The first generation TV Streamer is not compatible with: ReSound LINX™, ReSound LINX²™, ReSound ENZO™, ReSound ENZO²™, and ReSound Enya™ (ITE, ITC)



ReSound helps people rediscover hearing with solutions that emulate the natural function of the human ear. Our aim is to give your patients the opportunity to live a rich, active and fulfilling life unaffected by hearing loss.

We are the only hearing manufacturer that offers a complete ecosystem of integrated wireless connectivity across the entire product range, so your patients can stream audio, personalize their sound and control their hearing aids without the need for intermediary devices.

resoundpro.com/wirelessconnectivity

**RESOUND
NORTH AMERICA**

8001 Bloomington Freeway
Bloomington, MN 55420
1-800-248-4327
resoundpro.com

**RESOUND
CANADA**

303 Supertest Road
Toronto, Ontario M3J 2M4
1-888-737-6863
resoundpro.com

**RESOUND
GOVERNMENT SERVICES**

8001 Bloomington Freeway
Bloomington, MN 55420
1-800-392-9932
resound.com/governmentservices

© 2016 GN ReSound A/S. All rights reserved. ReSound is a trademark of GN ReSound A/S. Apple, the Apple logo, iPhone, iPad, iPod, iPod touch and Apple Watch are trademarks of Apple Inc., registered in the U.S. and other countries. App Store is a service mark of Apple Inc. Android, YouTube and Google Play are trademarks of Google Inc. The Android robot is reproduced or modified from work created and shared by Google and used according to terms described in the Creative Commons 3.0 Attribution License. Dolby and the double-D symbol are registered trademarks of Dolby Laboratories. Bluetooth is a trademark of Bluetooth SIG, Inc.

ReSound

rediscover hearing



MK604655 Rev B 2016.03