

Love the experience - introducing Discover

Bill Christman, AuD
Senior Regional Trainer

unitron Love the experience

A Sonosco brand

© 2019 Unitron. All rights reserved.

Learning objectives

After this course learners will be able to:

1. Identify the benefits of Unitron's newest platform and how the advanced performance can provide an easy first fit, sound performance patients love, and technology for busy lifestyles.
2. Identify the patient benefits of Unitron's first Lithium-ion rechargeable solution.
3. Describe the benefits of SoundCore technology on the new platform.



Discover



Fit

Wow from the get go



Discover – a better first fit experience



Better first sound

new initial frequency response



Faster fitting process

5-click workflow in TrueFit



More patient engagement

FLEX real-word assessment



Discover

Great style


Tailored fit

Stellar sound







Love

Ignite every conversation

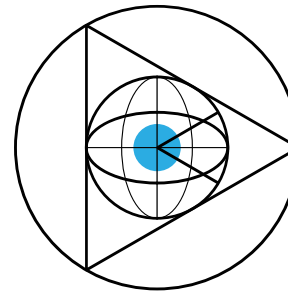


SoundCore



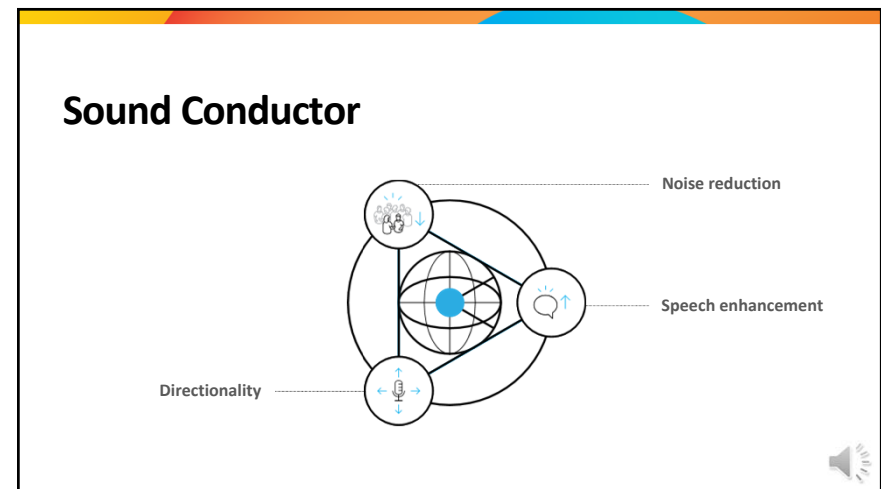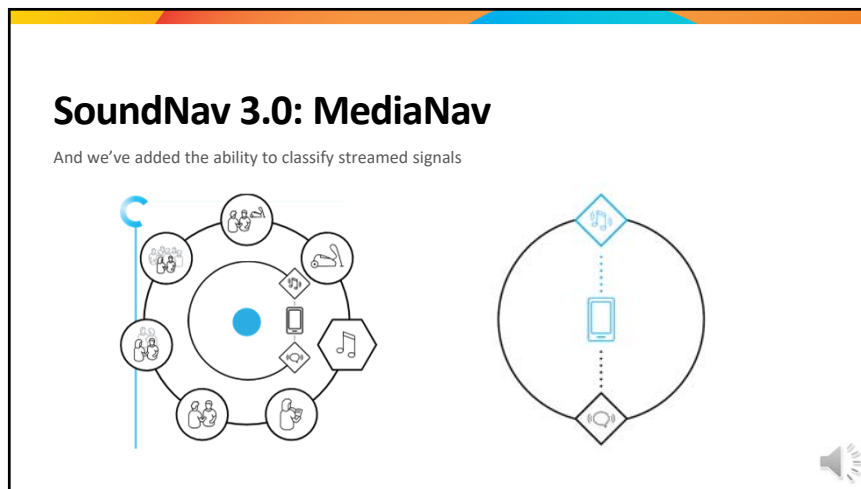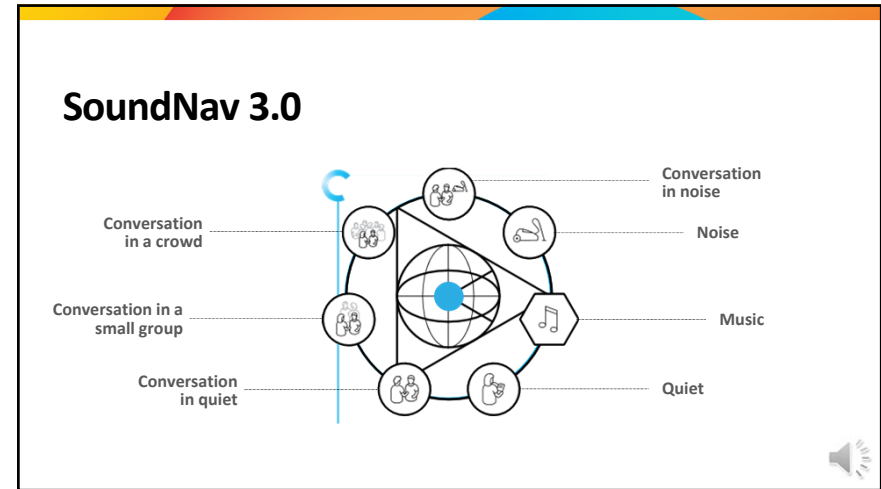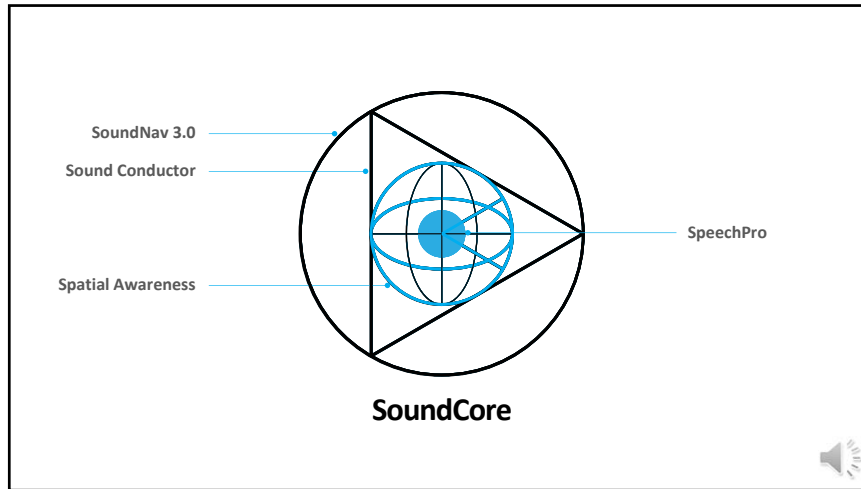


Directionality


Different directional strategies are applied depending on the classified environment

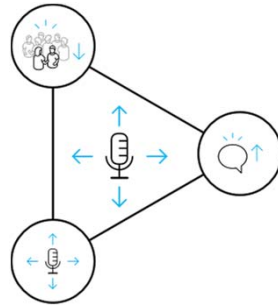


Discover - 4 technology levels

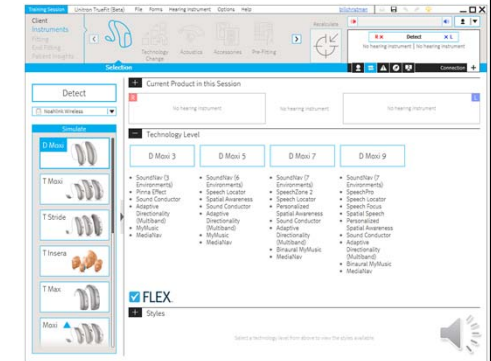


Discover 9

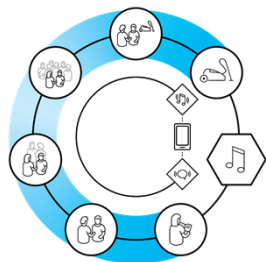


SoundNav 3.0

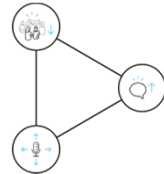

Sound Conductor

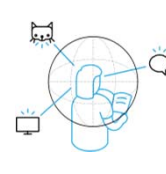


Spatial Awareness

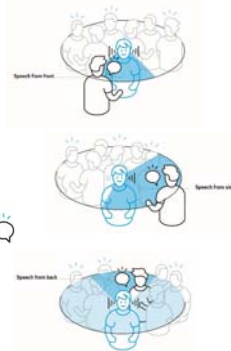


SpeechPro


Discover 7

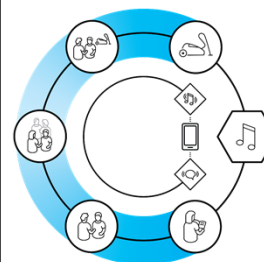


SoundNav 3.0

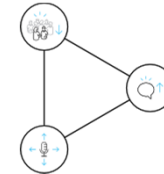

Sound Conductor

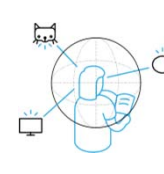


Spatial Awareness

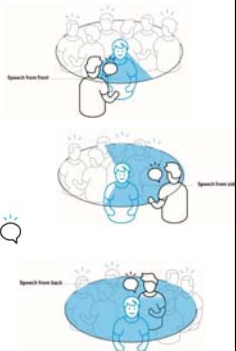


SpeechZone 2

Discover 5

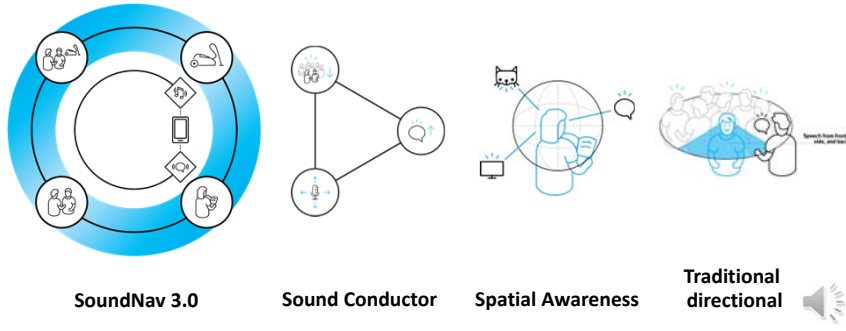

Discover 3

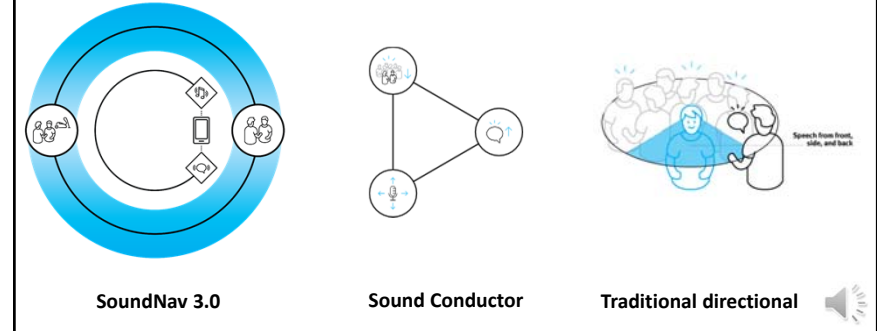

Make your patient feel 20 again

Discover

Ignite every conversation



Go

Chat. Charge. Repeat.



Fit.Love.Go.



Direct connectivity

Moxi Jump R



Moxi Fit

Brand new styles



Li-ion rechargeable



SWORD 3.0 chip supports direct iOS and Android connectivity



Bluetooth Classic
Hands-Free Profile (HFP) protocol
hands-free phone calls



Bluetooth® Low Energy (LE)
protocol for Remote Plus app,
fitting and Duolink



Bluetooth Classic
A2DP Profile protocol for direct,
binaural media streaming



AirStream Sonova propriety
protocol for audio streaming with
the TV Connector

*Smartphones and traditional mobile phones with a compatible Bluetooth Hands-Free Profile



Connectivity



ENTERTAINMENT:

Music
TV
Podcasts
YouTube
Netflix
Audiobooks

COMMUNICATION:

Phone calls
Skype
FaceTime
WhatsApp

INTERACTION:

Ratings
Virtual assistants like Siri
and OK Google



Impeccable styles





Moxi Jump R

- Automatically turns on when removed from the charger, turns off when inserted in the charger
- Turn on/off manually by holding the lower part of the button for 3 seconds
- Robust design with built-in, sealed design
- All day use on a full charge




Stock mode


To exit stock mode prior to fitting, place instruments in the charger and keep them inserted at least until the indicator light blinks after about 30 seconds




New receivers and domes

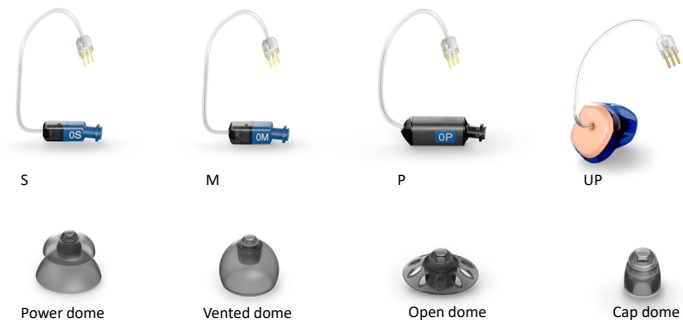

RIC Performance profile

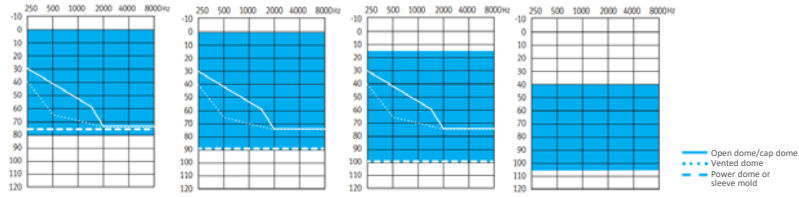

S receiver	M receiver	P receiver	UP receiver
111/46	114/50	122/58	130/67
Mild to moderate	Mild to severe	Mild to severe	Moderate to profound

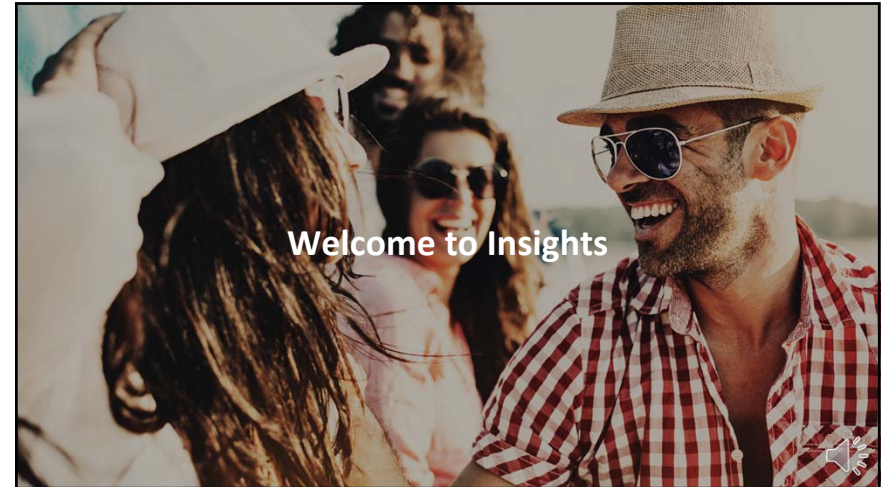

Insights



Log It All



Usage



Ratings



Overall satisfaction

unitron TrueFit.

unitron. myInsights



Say goodbye to the
status quo.
And get ready to love
the experience.

