

#### **Continuing 4S Foundation**

Sonic is dedicated to improving life through enhanced hearing by constantly focusing on our 4S Foundation

Sound that's natural

Speech understanding in noise

Simplicity in all we do

Style that stands out

With Captivate, that adds up to a better hearing experience

www.sonici.us



# **Revolutionary SoundDNA Platform**



From clarifying conversations to streaming online entertainment, Captivate immerses patients in the sounds of life.

SoundDNA Platform




#### **SoundDNA**

#### 2019 SoundDNA Updates

- More Processing Power
- New! Algorithm
- New! Setting
- New! Feature
- New! Benefits
- More Robust Offerings Overall





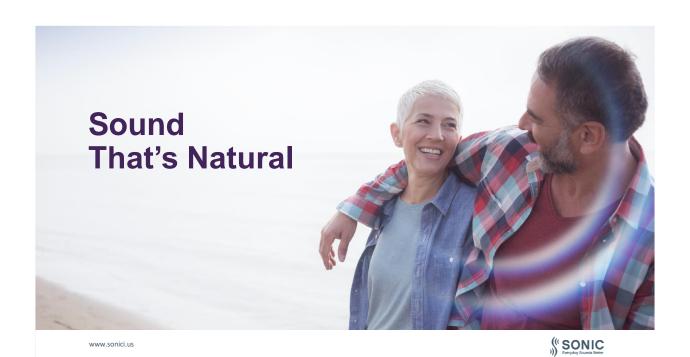

#### **Listen with Ease**

With Captivate, Sonic delivers a natural sound with several technologies in play

**SmartCompress** 

Adaptive Feedback Canceller Pro

Frequency Transfer

Extended Dynamic Range

**SmartMusic** 

www.sonici.us



# A More Open Fit With Less Risk of Feedback


**New!** Algorithm exclusive to 2019 Products on the SoundDNA Platform

1 Feedback Management system with 2 technologies




#### **AFC Pro Advantage**

- Fewer distractions and interruptions from the annoyance of feedback in both stable and changing conditions
- A better target match, for improved audibility of soft speech sounds
- More open fittings for a natural, own-voice sound

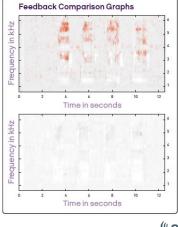

www.sonici.us



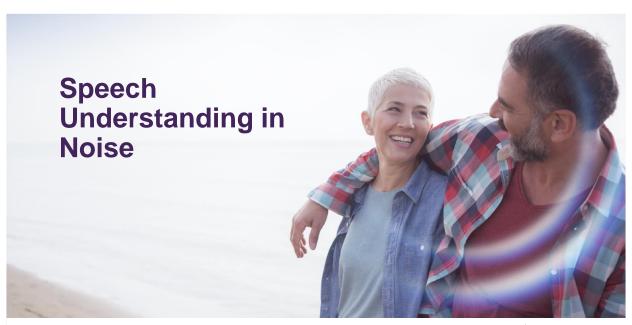

www.sonici.us

SONIC Everyday Sounds Better

#### **Hear Conversations More Clearly in Noise**

With Captivate, Sonic delivers what life sounds like in noise

**SPiN Management** 

Binaural Noise Management Impulse Noise Reduction Soft Noise Reduction Wind Noise Reduction

Tinnitus SoundSupport

www.sonici.us



#### **SPiN Management**

- SPiN Directionality
  - Uses null steering to select the microphones' polar patterns that will produce the optimal SNR in 16 independent frequency bands
- SPiN Noise Reduction
  - Designed to reduce the gain of noise compared to speech operating in 16 independent frequency bands
- SPiN Engage
  - New! Very High setting
    - Systems engage immediately at the onset of soft noise levels while the SNR is still high



#### **SPiN Management Advantage**


Better speech understanding in noisy listening environments

Settings that start at higher SNR levels to reduce background noise

Personalized to meet patients' preferences

www.sonici.us



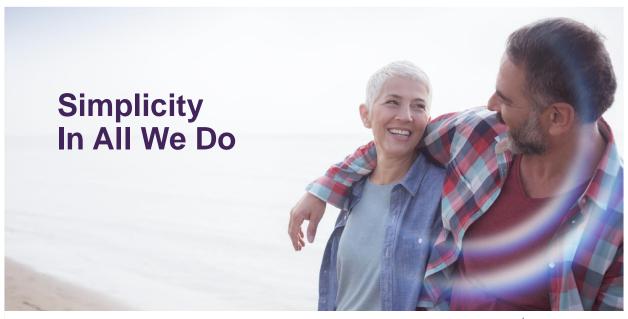

#### Make Life Easier with Captivate

Captivate delivers solutions that make sense

Lithium-ion Rechargeability

Dual-Radio System

If This Then That

EXPRESSfit® Pro

#### Recharge Your Day, Just Like That



New! Lithium-ion rechargeable battery solution

Robust battery life capacity for a full day's use following a short charging time



#### The Captivate Li-ion System

#### New model

miniRITE T R


Available in all Captivate levels
Charger is black - unbranded



www.sonici.us



#### Why Offer This Solution?

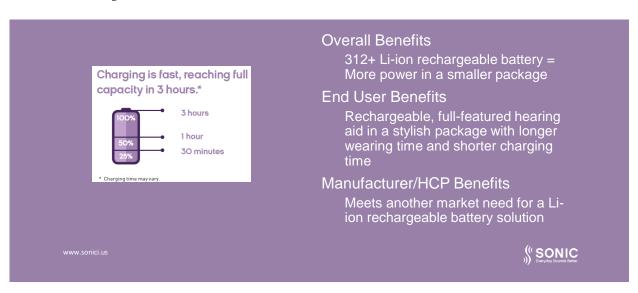


## EXPRESSfit® Pro Delivers...Again

Support for new Captivate products
Updated Transfer Fitting feature
Best Fit Fast SE improvements
Updated SPiN Management features
New! Battery Life Status feature





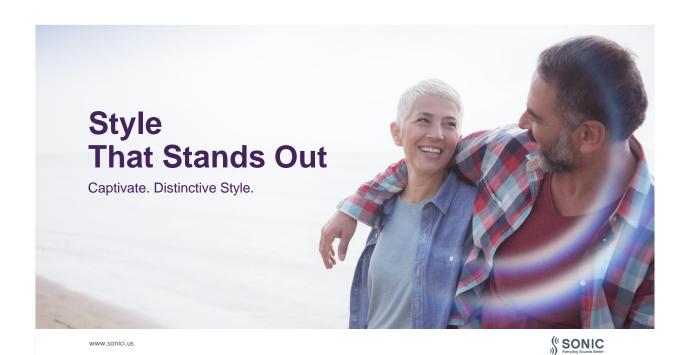

#### Flexibility in Focus

4 models miniRITE, miniRITE T, miniRITE T R, BTE 105

3 performance levels - 100|80|60

7 color options available



<sup>+</sup> Gold color available in the Captivate 100 miniRITE only.

### Design is in the **Details of Captivate**



- Dual Covered Microphones
   Allow for sophisticated directional features to enhance performance in
  - noise > Keep unwanted debris and moisture from entering the device
- Program Button
   Provides easy access for up to four listening memories
   Provides push-button volume control (miniRITE only)
- 3. Volume Control
- Programmable range
   Configurable alerting tones provide
   audible cue as volume is changed
- » NFMI enables binaural communication between devices » 2.4 GHz technology provides connectivity to external audio sources (mobile phone, television, computer, etc.)

- Battery Compartment: miniRITE, miniRITE T, BTE 105

- minifile 1, BIE 105

  » Easy-open design
  » Integrated on/off switch
  » Size 13 IBTE 105) or 312 [miniRITE,
  miniRITE I battery
  » Tamper-resistant lock keeps the
  battery secure from children

- Acoustic Options miniFIT System
  » Four receiver sizes in four lengths
  miniRITE\_miniRITE T, miniRITE T of the options
  (BTE 105 only)
  » Custom Molds and miniFIT Domes
  available in a variety of sizes and
  configurations
- 7. Rechargeability
- Nething and the solution in th

- Ingress Protection Rating

  » Hydrophobic coating and overall
  design protects against dust and
  continuous immersion in three feet of
  water or more
  » IP68
- ProWax miniFit System
   Protects receiver from moisture
   and cerumen (miniRITE, miniRITE T,
   miniRITE T R)
   Replaceable by the patient
- » Telecoil provides clear telephone reception (miniRITE T, miniRITE T R, BTE 105)
- 11. If This Then That (IFTTT) » Connects Captivate instruments to other internet devices via Sonic Go service
- - » Optional adapters for BTE 105 only

\* Not all features are available on all models—refer to the chart in the back of this guide for details









miniRITE

miniBITE T

BTE 10

#### Model Overview

				DIL 103
Battery Size	312	312	312	13
Power Model	Power Receiver	Power Receiver	Power Receiver	
Directional Microphones				•
Program Button	•		-	•
Volume Control	0	0	0	
Telecoil			•	•
Wireless Connectivity	•			
Receiver Sizes	Small (60) Medium (85) Power (100) Power (105)	Small (60) Medium (85) Power (100) Power (105)	Small (60) Medium (85) Power (100) Power (105)	
Earhook				
Thin Tube				
IP Rating	IP68	IP68	IP68	IP68
Tinnitus SoundSupport	•		-	•
FTTT Internet Connection*				
DAI/FM				
Li-ion Rechargeability				
ZPower Rechargeability	0			

■ Standard ■ Optional

<sup>\*</sup> Requires SoundLink 2 App.

